




Approved SleepTight Mouthpiece Cleaning Tablets

Fresh and Clean Mouthpiece Cleaning tablets can be purchased here.

[Fresh and Clean Dental Appliance Cleaning Tablets \(mouthguardstore.com\)](http://mouthguardstore.com)



Made in Germany

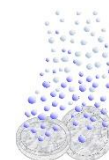

Description

- The Fresh and Clean dental appliance cleaning tablets will help keep your mouthguard or mouthpiece fresh and clean with a fresh minty taste.
- Stop the bacterial odor and stain build up on your oral device by using these safe, effective, and medically formulated sanitizing tablets.
- Clean your dental appliance at least once per week to keep your mouth healthy and feeling fresh!

Soft Dental Appliances such as SleepTight Mouthpiece Snoring appliances or sports mouthguards.

Instructions for Fresh and Clean

- Use a very soft toothbrush with a liquid antibacterial soap to remove any debris.
- Do not use toothpastes as they contain abrasives which can damage the surface over time causing it to be more susceptible to staining.
- Add the cleaning tablet to a glass of warm water (never hot water which could distort the device).
- Place the dental appliance into the effervescent solution and allow it to soak for 15 minutes.
- Remove and rinse and store dry. Keep the device away from pets as they may chew the device and make it unusable.